LIFEWAVE®

Y-age
glutathione

Y-age
carnosine

Y-age
aeon

Instructions

Vejledning | Anwendung | Mode d'emploi
Instruzioni | Instrucciones | Instrukcje
Pokyny | Arahan

INSTRUCTIONS FOR USE

Rotate through these point locations, using a Glutathione patch 5 to 7 days per week. You can use any of the point locations in this booklet for any Y-Age patch, but the preferred Glutathione locations are shown on page 7.

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Intended Use: To provide mild and temporary stimulation of acupressure points which improve the appearance of skin colour, texture and depth of wrinkles which combined with both Y-Age Aeon and Y-Age Carnosine. The device is a non-transdermal adhesive patch.

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INSTRUCTIONS FOR USE

Rotate through these point locations, using a Carnosine patch 1 to 3 days per week. You can use any of the point locations in this booklet for any Y-Age patch, but the preferred Carnosine locations are shown on page 13.

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Intended Use: To provide mild and temporary stimulation of acupuncture points which improve in energy flow, heart strength, flexibility, reduces stiffness and stress which supports a healthy inflammatory response. Also help to improve in the appearance of skin colour, texture and depth of wrinkles which combined with both Y-Age Glutathione and Y-Age Aeon. The device is a non-transdermal adhesive patch.

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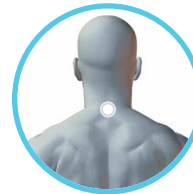
INSTRUCTIONS FOR USE

Rotate through these point locations, using an Aeon patch 5 to 7 days per week. You can use any of the point locations in this booklet for any Y-Age patch, but the preferred Aeon locations are shown on page 19.

When combining Y-Age products, use any 2 patches on any 2 of the Y-Age point locations in this booklet, applying the patches on either the midline or right side of the body.

Apply the patches to clean, dry skin in the morning. Patches may be worn for up to 12 hours before discarding. Keep well hydrated while using this product.

Intended Use: To provide mild and temporary stimulation of acupressure points which promotes relaxation, reduces stiffness and stress which supports a healthy inflammatory response, improve in the appearance of skin colour, texture and depth of wrinkles which combined with both Y-Age Glutathione and Y-Age Carnosine. To provide pain relief when Aeon is combined with a set of IceWave patches. The device is a non-transdermal adhesive patch.

Warnings: Remove immediately if you feel discomfort or skin irritation occurs. Do not reuse patch once removed from the skin. For external use only. Do not ingest. Do not use on wounds or damaged skin. Ask a health professional before using if you have a health condition, any questions or concerns about your health. Do not use if pregnant or nursing. Product to be stored at room temperature.

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LifeWave Europe Ltd., Raheen Ind. Est., Athenry, Co. Galway H65 NW84, IRL.
LifeWave, Inc., 9775 Businesspark Ave, San Diego, CA 92131, USA.